

Accessories to meet your individual requirements



For higher hygiene demands a paper roll holder can be mounted easily at the chair's back segment.




The practical reading lamp offers additional light for your patients, e.g. for reading.



The IV-pole can be mounted space-saving on the left or right side of the base frame.



If desired the SlimLine can be equipped with an equipotential bonding.

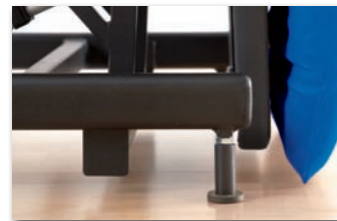

If necessary the stable feet allow an even lower entrance height of only 54 cm.



The practical bag holder can be mounted to the armrest fast and easily.



Higher flexibility and comfort for your patients offers the space-saving and easy to adjust table.



Available colours: Armrests, base upholstery and neck pillow are always effected in Anthrazit. Standard upholstery colour is Atoll. Many other colours are available. Colourings by jeans or other textiles are excluded from any manufacturer's guarantee. Please ask for our colour chart.

Technical data:

Dimensions

Seat segment width	54 cm
Chair width, incl. armrests	80 cm
Height (fixed, with wheels, standard)	57,5 cm
Height (fixed, with feet, optional)	54 cm
Total length	185 cm
Total weight	51 kg

Safety

Safe maximum load	175 kg
Angle of head down position	- 10°

Adjustments

Back segment angle	37° – 90°
Seat segment angle	- 4° – 20°
Leg segment angle	0° – 90°

Electrical specifications

Protection class	2
Protection category	IPX 6
Protection category (motors)	IPX 4
In compliance with	EN 60601
Motors	2



Bionic Medizintechnik GmbH ... a JMS-corporation · Max-Planck-Straße 21 · D-61381 Friedrichsdorf/Germany
Phone: +49(0)6172-75760 · Fax: +49(0)6172-757610 · info@bionic-jms.com · www.bionic-jms.com

Right for technical modifications reserved
B-MOBSLE/04.2017-02

Space-saving, flexible and comfortable - SlimLine for the infusion therapy



TekMed
1300 720 727
www.tekmed.com.au

TekMed NZ
0800 456 217
www.tekmed.co.nz



FS521104



SlimLine – minimal space requirements – maximum comfort and flexibility for your patients



Easy entrance and exit from the chair's front is possible due to the low seat height of only 57,5 cm, the adjustable right angled leg segment and the perfect seat inclination.



The SlimLine Therapy Chair has been especially developed for the infusion therapies and can be adjusted with two motors into a convenient comfort position fast and easily.



In case of an emergency the SlimLine can be moved into a lowered head position for cardio pulmonar rescue by pushing a button. Although the base frame is small and space-saving it offers perfect stability.

The SlimLine Therapy Chair has been developed for the needs in your oncology clinic and the infusion therapies.

Also in small rooms the slim and light Therapy Chair offers with its two motors the necessary flexibility and highest comfort to guarantee your patients a pleasant treatment.

With the help of the user-friendly and clear to understand remote control the back, seat and leg segment can be adjusted into the desired position. Here the leg segment moves synchronously with the seat, the back segment can be adjusted separately. By pushing a button the SlimLine is moved into a relaxed comfort position or in case of an emergency into a lowered head position. Due to the practical holder on each armrest the remote control is always accessible easily by the patients.

The special concept of the armrest guarantees highest stability also during entrance and exit from the chair's front. The armrest can be lowered fast and easily by pushing a button and therefore fit the patient's position.

Special comfort is guaranteed by the removable and easy to clean mattress in combination with the unique base upholstery. The ergonomically shaped back segment and the separately upholstered seat allow a comfortable and relaxed position during long term treatments.

The comfortable neck pillow can be adjusted easily into the desired position with a counterbalance.

Additional individual accessories such as IV-pole, paper roll holder or table can be mounted space-saving directly on the SlimLine. If your patient wishes to read or work during treatment, there is an optional mountable reading lamp. Equipotential bonding can be ordered easily.

To fit to your environment the colour of the SlimLine mattress can be chosen from more than 30 different colours. The base frame, the armrests and the optional neck pillow are always effected in Anthrazit.



The neck pillow offers higher comfort and can be adjusted individually with a counterbalance.



More flexibility: The SlimLine is equipped with wheels (7.5 cm ø) with single brake.



The illustrated remote control has two pre-programmed buttons for a simple operation. Each armrest is equipped with an easy to reach holder for the remote control.



The stable armrests can be lowered by pushing a button and therefore adjusted to the needs of your patients.



The unique base upholstery consists of an ergonomically shaped back segment and a flexible seat segment.



The removable and easy to clean mattress offers in combination with the base upholstery long-lasting comfort.